



LMPerioKit by Style Italiano

LM-PerioKit by Style Italiano is a set of periodontal instruments chosen by Dr. Giano Ricci.

LM periodontal instruments are ultimate in metallurgy. The high wear resistance and hardness of the blade metal makes them efficient, long-lasting and comfortable. The ergonomic design with a LM-ErgoSense silicone handle gives the instruments a non-slip and natural grip. Good tactile sensitivity and easy handling of the instruments guarantee the best clinical results. The LM-ErgoSense handle found to have the highest rated ergonomics in studies by the Finnish Institute of Occupational Health.

LM-ErgoSense handle

- Great tactile sensitivity
- Larger diameter & optimal design
- Unique intelligence with RFID technology
- Rated best in clinical testing



LM-PerioKit LM 6842 ES

Includes the following instruments:

North Carolina Probe LM 54B ES

For detecting calculus, cavities and abnormalities on the tooth surfaces as well as probing the depth of pockets. With 5-10-15 mm markings.

Furcation Probe LM 20B-21B ES

For examination of furcation lesions from different angles in both the upper and lower jaws.

McCall 17S-18S LM 279-280 ES

For the removal of sub- and supragingival calculus. Excellent in tight interdental spaces.

Concavity Diamond File SV 5-6 LM 255-256D ES

For acute furcation lesions. Diamond coated round curved toe.

Mini Gracey 7/8 LM 207-208M ES

For premolars and molars. For buccal, mesial and lingual surfaces. Great for deep and narrow pockets.

Mini Gracey 11/12 LM 211-212M ES

For premolars' and molars' mesial surfaces. Great for deep and narrow pockets.


Mini Gracey 13/14 LM 213-214M ES

For premolars' and molars' distal surfaces. Great for deep and narrow pockets.

Implant Mini Universal curette LM 283-284MTI EM

Titanium tips for sensitive but still effective implant maintenance. Universal design works well in tight interdental spaces. With replaceable LM-ErgoMix tips.

“The LM-PerioKit by Style Italiano contains the essentials for the best possible therapy”

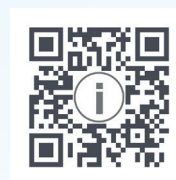

For more details about LM ErgoSense handle, scan the QR and visit our website.

Instruments have been selected in the spirit of Style Italiano whose philosophy is to make treatment feasible, teachable, repeatable. Periodontal therapy too can follow this approach. It's my feeling that, in the best interest of his patients, every practitioner should be able to properly execute non surgical therapy, crown lengthening procedures, free gingival grafts, simpler surgical root coverages.”

Giano Ricci
MD. DDS. MScD.

